ÓRTESE CORRECTIVA DUPLA DO HALLUX VALGUS

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DOUBLE PROTECTION HALLUX VALGUS CORRECTIVE ORTHOSIS

INDICATIONS

Daily use of the EPITACT® double orthosis helps correct hallux valgus deformity and relieve associated pain. In addition, it acts effectively against plantar pain.

CONTRA-INDICATIONS

- Big toe that does not move sideways.
- Do not use on unprotected wound.
- Persons suffering from arthritis, diabetes, polyneuritis, lower limb neuropathy, varicose ulcers or severe oedema: never use the product on a wound, monitor the condition of your foot during use, and in all cases seek advice from your doctor or chiropodist.

PRECAUTIONS FOR USE

- Do not use this product during the night.
- There are two versions: left foot and right foot. The black part of the device should be in contact with the insole.
- Wash before first use.

- PLACEMENT:

- 1 Slip the strap onto the foot.
- 2 Place the cap on the big toe and align the marker with the bunion. Insert your second or third toe into the elastic of the plantar protection according to your preference.
- 3 If the toe is still misaligned, gradually tighten the strap until the toe is straightened.

During the first few days of use, you can apply a light correction and gradually increase it to get used to the device.

CLEANING


- Washing after use at 30°C in the washing bag supplied or wash by hand (without rubbing).
- Discard and replace the device as soon as cracks in the gel or fabric or seam damages occur.
- Leave to air-dry. Do not iron.

COMPOSITION: 54% silicone, 46% polyamide/elastane fabric.

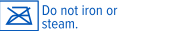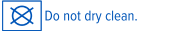
In case of pain or any other unusual reaction, contact your doctor. The silicone tendon in contact with the skin may cause mechanical irritations, such as blisters (< 0.1%). Any major incidents involving the device must be reported to the manufacturer and the relevant health authority. MILLET Innovation cannot be held responsible for failure to use the product in accordance with the manufacturer's recommendations. French is the reference language. CE initial marking: 2012.

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